

More Accurate Methods for Detecting and Quantifying Fat from Magnetic Resonance Images

View U.S. Patent No. 8,000,769 in PDF format.

WARF: P07150US

Inventors: Scott Reeder, Huanzhou Yu

The Wisconsin Alumni Research Foundation (WARF) is seeking commercial partners interested in developing methods for improving the accuracy of fat measurements from MR images by removing the effects of noise during the quantification process.

Overview

Quantifying the amount of fat in the liver is crucial for the detection of non-alcoholic fatty liver diseases including steatosis, fibrosis, cirrhosis and liver failure. Accurate fat measurements can lead to early treatment and prevention of liver diseases. Currently, the best method for the measurement of liver fat is a biopsy, in which bits of the liver are surgically removed and analyzed.

MRI techniques have been developed to image fat, but currently are only useful for fat detection. Because these techniques produce large amounts of noise, the quantification of fat cannot be accomplished reliably. An accurate, non-invasive fat quantification method would be desirable.

The Invention

UW-Madison researchers have developed two methods for improving the accuracy of fat measurements from MR images. These methods remove the effects of noise during the quantification process. Large amounts of noise leads to an over- or underestimate of the amount of fat. The removal of this noise results in a more accurate fat image.

The two methods are known as the phase constrained and signal normalization methods. The phase constrained method assumes water and fat to have a common, equal phase. This method calculates the fat fraction directly from real water and fat components, as opposed to complex signal components, to remove noise bias. The signal normalization method takes an image signal when the water and fat are in phase. Two or more other image signals are obtained and normalized using the previous in phase image signal to remove any noise bias.


Applications

- Measuring and quantifying fat from MR images
- · Detecting fat in the liver, skeletal muscles, adrenal glands and other tissues or organs


Key Benefits

Improves the accuracy of quantifying fat from MR images

We use cookies on this site to enhance your experience and improve our marketing efforts. By continuing to browse without changing your browser settings to block or delete Additional information cookies, you agree to the storing of cookies and related technologies on your device. See our privacy policy

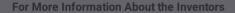



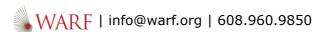

Scott Reeder

Related Intellectual Property

• View Divisional Patent in PDF format.

Tech Fields

• Medical Imaging: MRI

For current licensing status, please contact Jeanine Burmania at jeanine@warf.org or 608-960-9846

