



## Robust Biological and Chemical Detection Method and Microfluidic Device with Liquid Crystal Sensing Element

[View U.S. Patent No. 8,925,373 in PDF format.](#)

**WARF: P09136US**

Inventors: Hongrui Jiang, Daming Cheng

**The Wisconsin Alumni Research Foundation (WARF) is seeking commercial partners interested in developing an improved biological and chemical detection device and method integrating liquid crystals.**

### Overview

Detecting chemical and biological agents at low levels is important in military applications, drug enforcement, airport security and environmental monitoring. A variety of devices and methods currently are employed for this purpose. For example, gas chromatography has been implemented widely into airport security since 2001 for chemical detection of explosives, and nucleic acid detection methods have been employed by U.N. officials to identify biological weapons. However, most current detector technologies require advanced technical training and are impracticable for use in the field due to large size or fragility.

UW-Madison researchers previously developed a simple, robust device and method for detecting biological agents in water or similar aqueous substances (see WARF reference number P06457US). The device employs a thin polymeric membrane that dissolves in the presence of a target compound at a speed correlated to the concentration of the target.

Liquid crystals also have been used for the detection of biological or chemical agents. Liquid crystals exhibit unique anisotropic physical properties that make them particularly well suited for detecting trace amounts of compounds.

### The Invention

UW-Madison researchers have now developed an improved method for autonomously generating stabilized liquid crystal thin films and two microfluidic devices that employ the technique for detecting trace amounts of biological agents and chemical compounds. The handheld microfluidic devices each contain a microchannel defined by grooved polymer materials sandwiched between glass substrates. Priming the device involves filling the microchannel with liquid crystal material, which fills specific nickel-plated structures in the channel, and flushing the liquid crystals outside the container with the laminar flow of an aqueous solution. This method allows for automatic formation and rapid regeneration of the stable aqueous/liquid crystal interface.

In the presence of a target compound the orientation of the liquid crystals changes, altering optical properties of the liquid crystals through a phenomenon known as optical birefringence. After the analyte has been introduced into the channel, a white light is passed through a first polarizing lens, the microfluidic device, and a second orthogonally oriented polarized lens. The intensity of the light, determined by the degree of optical birefringence, is detected by a microscope to confirm the presence or absence of the specific target in the aqueous solution.

We use cookies on this site to enhance your experience and improve our marketing efforts. By continuing to browse without changing your browser settings to block or delete cookies, you agree to the storing of cookies and related technologies on your device. [See our privacy policy.](#)



advanced technical training, the improved detection apparatus will greatly enhance in-field applicability of biological and chemical sensor technologies.

## Applications

- Biological and chemical sensing at very low concentrations
- Early detection of biological and chemical agents in military applications
- Industrial chemical hazard warning systems
- Drug enforcement
- Airport and border security
- Environmental pollution monitoring

## Key Benefits

- Enables automatic operation of integrated liquid crystal sensing devices
- Eliminates need for careful manual handling and technical training
- Provides real-time, repeatable sensing tasks in the field
- Utilizes liquid crystal sensing within the microfluidic device
- Allows automatic formation of sensing interfaces
- Provides control over the target phase/liquid crystal interaction

## Stage of Development

Prototype of both devices were able to detect adsorption of the cationic surfactant dodecyltrimethylammonium bromide (DTAB) as well as biomolecular interactions between a phospholipid and the enzyme phospholipase A<sub>2</sub> (PLA<sub>2</sub>).

## Additional Information

### For More Information About the Inventors

- [Hongrui Jiang](#)

### Related Technologies


- [For more information about detecting and determining the concentration of a biological agent in liquids, see WARF reference number P06457US.](#)

### Tech Fields

- [Analytical Instrumentation, Methods & Materials : Biodefense](#)
- [Analytical Instrumentation, Methods & Materials : Microfluidics](#)
- [Clean Technology : Monitoring, remediation & waste reduction](#)

For current licensing status, please contact Jennifer Gottwald at [jennifer@warf.org](mailto:jennifer@warf.org) or 608-960-9854

## Figures

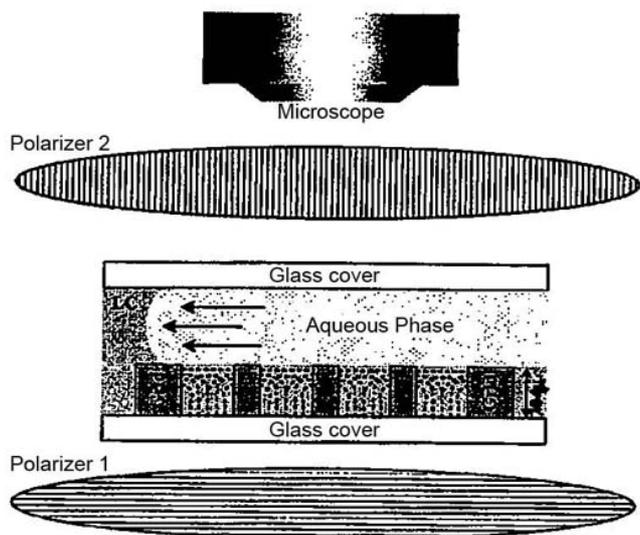

We use cookies on this site to enhance your experience and improve our marketing efforts. By continuing to browse without changing your browser settings to block or delete cookies, you agree to the storing of cookies and related technologies on your device. [See our privacy policy.](#)

OK



WARF  
Wisconsin Alumni Research Foundation

| [info@warf.org](mailto:info@warf.org) | 608.960.9850

We use cookies on this site to enhance your experience and improve our marketing efforts. By continuing to browse without changing your browser settings to block or delete cookies, you agree to the storing of cookies and related technologies on your device. [See our privacy policy.](#)

OK



**WARF**  
Wisconsin Alumni Research Foundation

| [info@warf.org](mailto:info@warf.org) | 608.960.9850